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TABLE OF CONTENTS

SAFETY PRECAUTIONS	26
PRECAUTIONS BEFORE USE	29
PRECAUTIONS IN USE	31
MAINTENANCE AND INSPECTION PRECAUTIONS	32
CONTENTS OF PACKAGE	33
NAME OF PARTS	35
HOW TO CHARGE	36
EXPLANATION OF FUNCTIONS	37
HOW TO USE	39
Q & A	44
TROUBLESHOOTING	45
PRODUCT SPECIFICATIONS	46
RECYCLING OF LITHIUM-ION BATTERY	47
ASSURANCE PROVISION	48
WARRANTY CARD	

SAFETY PRECAUTIONS

Be sure to read this instruction manual carefully before using the product so that you are able to use it properly at all times.

To prevent injury, loss of life and damage to property, the following symbols are used.

 DANGER	Indicates an imminently hazardous situation which, if not avoided, could result in death or serious injury.
 WARNING	Indicates a potentially hazardous situation which, if not avoided, could result in death or serious injury.
 CAUTION	Indicates a potentially hazardous situation which, if not avoided, could result in minor injury or property damage*. * Property damage means damage to house, household belongings, etc.

■ MEANINGS OF SYMBOLS

	This symbol indicates a prohibition. * The details appear inside or near the symbol.
	This symbol indicates a requirement that must be followed.
	This symbol indicates a caution.

SAFETY PRECAUTIONS

To avoid accidents and discomfort

DANGER

 <p>Prohibition</p>	<ul style="list-style-type: none"> ◆ Persons using the following medical electronic equipment must not use this product. <ul style="list-style-type: none"> • Malfunction of medical electronic equipment may cause injury or discomfort. • Biological implant such as pacemaker • Life support electronic equipment such as artificial heart-lung machine • Wearing-type medical electronic equipment such as electrocardiograph or infuser
------------------------------------------------------------------------------------------------------	------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------

DANGER

 <p>Prohibition</p>	<ul style="list-style-type: none"> ◆ The following persons must not use this product. <ul style="list-style-type: none"> • Persons with heart disease, abnormal blood pressure or malignant tumor • Persons with a febrile disease, infectious disease or acute disease • Persons with a weakened immune system or who use immune suppressants due to immunity disorders such as AIDS • Persons with hemophilia, a bleeding disorder, diabetes or facial neuralgia • Persons with paresthesia or dysalgia • Persons with Keloid diathesis, facial melanosis, photosensitivity or a pigmentation disorder • Persons with capillary dilatation due to long-term steroidal hormone treatments or liver dysfunction • Persons who have allergies, atopic dermatitis, skin diseases • Persons who experience epileptic seizures or discomfort due to flashing lights • Persons who have aftereffects of head injury or ophthalmologic affects due to traffic accidents • Persons who are pregnant, menstruating or breastfeeding • Infants, growing children and persons who cannot express their feelings in words • Persons who are under the influence of alcohol or experiencing a hangover • Persons who take drugs that cause drowsiness • Persons who are extremely tired or do not feel well • Persons who have had skin problems such as spots or rashes due to artificial fibers • Persons who have or have had skin inflammation due to cosmetics • Persons who have had chemical peeling treatment at a beauty salon, etc. in the immediately recent past • Persons who have previously had skin problems after using beauty equipment
 <p>Prohibition</p>	<ul style="list-style-type: none"> ◆ Do not use the product on the following body parts or areas. <ul style="list-style-type: none"> • Except specified body parts • Around heart, eyeballs, upper eyelids, lips and mucous membranes such as inside the mouth • Body parts with morbid spots, moles and viral warts • Body parts with wounds, eczema, psoriasis, herpes, purulent pimples • Body parts with pathological spots due to itching, flushing or physical stimulus

SAFETY PRECAUTIONS

To avoid accidents and discomfort

WARNING

Compulsory

- ◆ The following persons must consult a doctor before use.
 - Persons who are seeing a doctor
 - Persons who have had serious diseases
- ◆ If you feel any pain or abnormality on or in your skin or body, stop using the product immediately
- ◆ Please follow the instructions for use. Otherwise, the product in operation may cause adverse effect on health.
- ◆ Use the product within the recommended period of time (30 minutes per day), unless otherwise permitted. Application on the same region of skin for longer than the specified time may cause potential damage to skin.
- ◆ Do NOT avoid or deactivate safety mechanisms specific to this product.

PRECAUTIONS BEFORE USE

WARNING

Main Body

 Prohibition	<ul style="list-style-type: none">• Do not use this product at the same time as using other beauty equipment.• Do not use this product, which is specialized for use as a beauty product, for any other purposes, such as therapeutic purposes. Do not use this product under conditions of breakdown or failure.• Do not look directly or indirectly at the LED light.• Do not use on pets.• Do not use near flammable substances such as thinner, benzene, polish removers or sprays.
 Prohibition of use in wet areas	<ul style="list-style-type: none">• Do not use or keep in the bathroom or high-humidity areas. Do not leave anywhere where contact with water can readily occur, such as on a sink.
 Prohibition against getting wet	<ul style="list-style-type: none">• Do not immerse in water, pour or spill water over or wash with water (risk of electric shock, fire or failure).
 Compulsory	<ul style="list-style-type: none">• Before using the product, cleanse your skin and apply an appropriate moisturizer.

PRECAUTIONS BEFORE USE

WARNING

POWER

 Prohibition	<ul style="list-style-type: none">• Do not use with any AC adapter and USB cable other than those provided with the product. Do not connect the dedicated adapter and USB cable to any other equipment (risk of fire, electric shock or failure).• Be sure that the USB port is not contaminated (risk of fire, electric shock or failure).• Do not damage the USB cable.* Do not damage, bend more than necessary, pull, twist, place heavy objects on top of or interpose (risk of fire, electric shock or damage).• Only use when firmly plugged into a power outlet.• Do not plug into an overloaded power outlet (risk of heat generation or fire).
 Prohibition of use in wet areas	<ul style="list-style-type: none">• Do not use or keep in the bathroom or high-humidity areas. Do not leave anywhere where contact with water can readily occur, such as on a sink.
 Prohibition against getting wet	<ul style="list-style-type: none">• Do not immerse in water, pour or spill water over or wash with water (risk of electric shock, fire or failure).
 Prohibition of use with wet hands	<ul style="list-style-type: none">• Do not plug in or unplug the AC adapter with wet hands (risk of electric shock or injury).
 Unplug	<ul style="list-style-type: none">• Be sure to hold onto the AC adapter, not the USB cable, when plugging in or unplugging.• Do not use if the AC adapter is damaged or not firmly plugged into a power outlet.
 Compulsory	<ul style="list-style-type: none">• Be sure to use only household outlets.• Make sure that the AC adapter and the USB cable are fully and firmly plugged in to the outlet or port.

PRECAUTIONS IN USE

WARNING

 Unplug	<p>In the event of the following, stop using the product, turn off the power and unplug the AC adapter.</p> <ul style="list-style-type: none">• Earthquake or lightning• Smoke, abnormal heat, odor or distortion coming from or affecting the main body, AC adapter or the USB cable <p>* Immediately consult your dealer or our customer support center.</p>
 Prohibition	<ul style="list-style-type: none">• During use, this product may become slightly warm. This is not a safety issue.• Do not bring this product into contact with metal objects or any other electrical appliances.
 Compulsory	<ul style="list-style-type: none">• Charge to full capacity before using for the first time.* When not sufficiently charged before use, the product may switch off during use.• Make sure to remove all metal objects and glasses during use.* Metal objects: earrings, necklaces, watches, etc.• If you feel any pain or abnormality on or in your skin or body during use, immediately stop using and consult a doctor.
 Compulsory	<p>Turn the product off and discontinue use in the following cases.</p> <ul style="list-style-type: none">• Smoke, abnormal heat, odor or distortion coming from or affecting the main body or the USB cable• Abnormal heat, swelling, or crack of the remote controller <p>* Immediately consult your dealer or our customer support center.</p>

CAUTION

Precautions on use of product Read the following carefully before use.

Skin condition depends on body condition.

Monitor the condition of your skin while using this product. If abnormalities such as rash, redness or itchiness appear on your skin, discontinue use and consult a dermatologist.

MAINTENANCE AND INSPECTION PRECAUTIONS

WARNING

 <p>Prohibition of disassembly</p>	<ul style="list-style-type: none">• Do not disassemble, modify or repair. For repair, consult your dealer or our customer support center.
 <p>Prohibition</p>	<ul style="list-style-type: none">• Keep the AC adapter and plug pins free of dust (risk of fire, electric shock or failure).* Clean the product regularly to prevent dust adhesion.• Do not drop or hit this product.• Do not heat or throw into a fire.• Keep out of reach of children and do not let children use this product (risk of accident or injury).• Do not store this product in the following places:<ul style="list-style-type: none">Places subject to direct sunlight, high humidity or dustPlaces subject to high or freezing temperatures or inside a car on a hot dayPlaces subject to risk of ignition, rain or windPlaces subject to chemicals, gas generation or excessive vibrationClose to heating appliances, appliances generating heat and sources of magnetic fields (risk of heat generation, fire or rupture due to the built-in rechargeable battery).
 <p>Compulsory</p>	<ul style="list-style-type: none">• How to clean<ul style="list-style-type: none">Wipe this product using a soft cloth that has been dampened with water. Wiping with a hard cloth may leave scratches on the surface or cause a decline in performance.* Do not use alkaline agents, alcohol, thinner, benzene or other solvents, polish removers, polishing powder, hot water or chemicals.* After every use, remove any residual beauty essence or gel remaining on the product.


INSPECTION

If this product remains out of use for a period of 6 months or more, be sure to inspect it for any abnormalities (e.g. LED lighting and cracks) before using it again.
Using under faulty conditions may result in an accident.

CONTENTS OF PACKAGE

Package 2: Mask Only

<Mask>

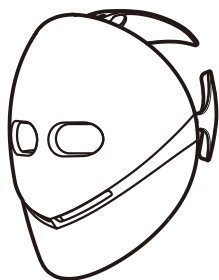
<Remote Controller
(with built-in battery)>

<Instruction Manual
with Warranty>


<USB Cable (long)>
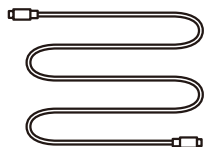
<AC Adapter>

<Dust Bag>

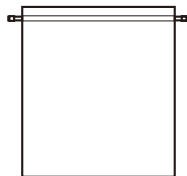Package 3: Neck Only


<Neck>


<Remote Controller
(with built-in battery)>

<Instruction Manual
with Warranty>
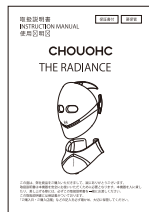
<USB Cable (long)>

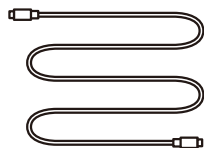
<AC Adapter>


NAME OF PARTS


<Mask (front)>

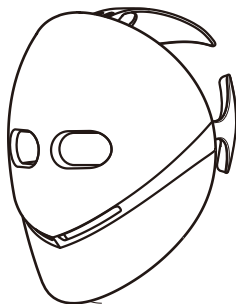
<Mask (back)>

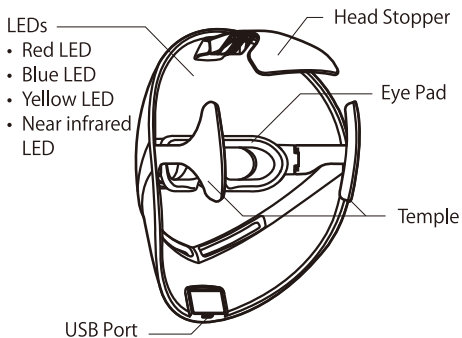
<Neck (front)>

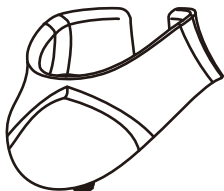
<Neck (back)>
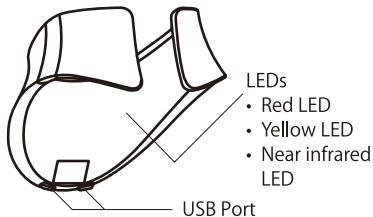
<Remote Controller>

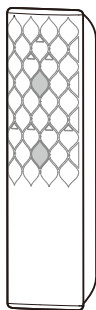USB Port

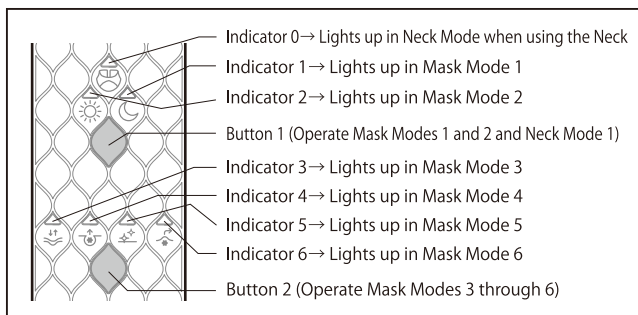
<AC Adapter>

USB Port

HOW TO CHARGE

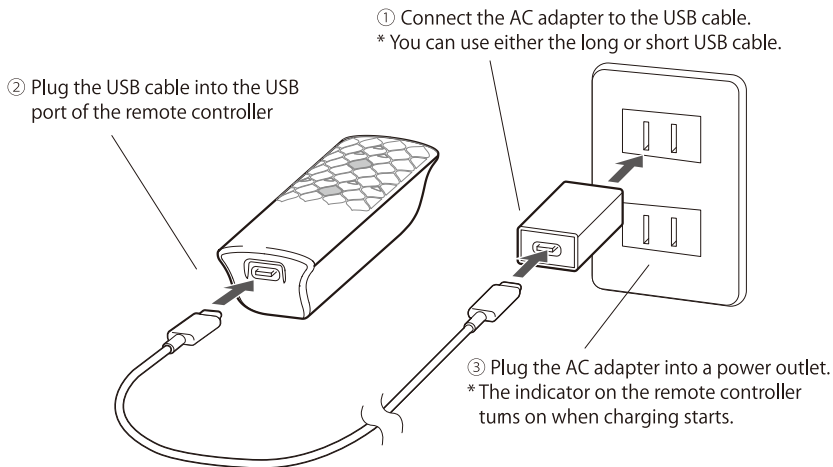- Charge to full capacity when the remaining battery life is low or before using the product for the first time.
- * The battery will gradually degrade with repeated charging and use.
The number of uses per charge may become less than indicated. This does not constitute failure.
- * The battery is not replaceable.
- * If you will not use the product for a long time, conduct a test operation once a month and charge to full capacity before storing.
- * When the battery goes below 10%, the indicator lamp begins to blink.

④ After charging completes, unplug the AC adapter from the outlet.

<Indicator>

Remaining battery: low	Indicator 0, Indicator 1 and Indicator 2 blink
While charging	<p>Indicator 0, Indicator 1, and Indicator 2 will light up in succession</p> <div style="display: flex; justify-content: space-around; align-items: center;"> <div style="text-align: center;"> <p>Indicator 0</p>  </div> <div style="text-align: center;"> <p>Indicator 1</p>  </div> <div style="text-align: center;"> <p>Indicator 2</p>  </div> </div>
Fully charged	Lights-out

Mode functions

Red LED and NIR LED reach both top and deep layers of the dermis, where they encourage collagen rebirth for the improvement of wrinkles and sagging in the neck.

Performing its functions all day long can make skin tired, and compromise its metabolism, accelerating skin aging unless remedied. This mode tightens the skin texture and pores, and helps improve sagging, an essential skincare step before bedtime.

In addition to germ-killing benefits for the skin surface and pores, and skin brightening benefits, this mode helps make wrinkles less noticeable by tightening the skin, achieved by stimulating fibroblasts to encourage collagen rebirth. The mode helps calm bacterial floras that have proliferated overnight, and awakens and tightens the skin for better makeup adhesion.

This mode is designed for greater skin tightness, achieved by encouraging collagen rebirth in both top and deep layers of the dermis. It is the mode to choose for addressing wrinkles and sagging.

Enlarged pores are due to overactive sebum glands, plus clogging caused by sebum plugs (sebum and dead skin cells). Pores can be tightened by controlling the activity of sebum glands attached to hair follicles, and by encouraging collagen rebirth in open pores.

Actions of yellow LED and red LED encourage epidermal cell metabolism leading to smoother removal from the epidermis of melanin buildup, thereby improving dark spots and dullness.

Causes of acne are sebum, sebum plugs, and proliferation of germs such as acne bacteria and staphylococcus bacteria. This mode controls acne activity by controlling bacterial proliferation, and by calming acne inflammation.

* The white circles in the "Duration" column of the table on the left signify that the corresponding LEDs are lit.

HOW TO USE


(Using the Mask and the Neck at the same time)

Please make sure this product is fully charged before use. Do not wear your glasses when you use this product.


* You cannot use commercially available USB cables. Use the USB cable provided with the product.

* When the Neck is connected to the Mask, the Neck will work in conjunction with the Mask.

- ① Plug the USB cable (long) into the USB port of the remote controller

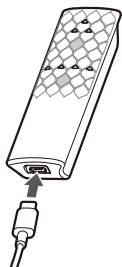
- ② Plug the USB cable (long) into the USB port of the Neck.

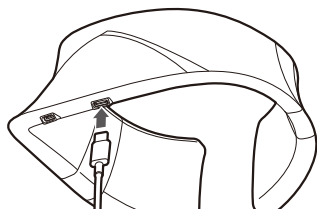
* You can plug the cable into either of the two USB ports.


- ③ Plug the USB cable (short) into the USB port of the Neck.

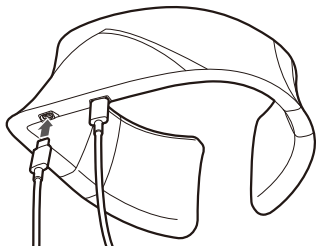
- ④ Plug the USB cable (short) into the USB port of the Mask.

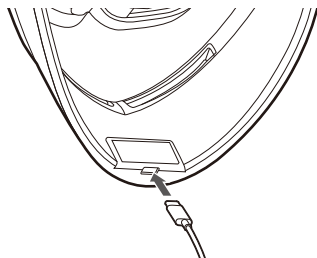- ⑤ Open the Neck with your hands to put it around your neck.

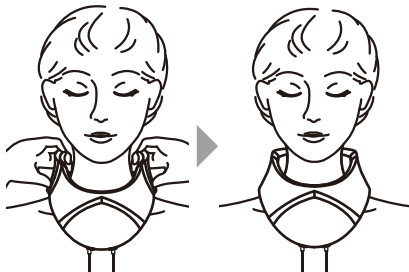
- ⑥ Lift up the head stopper. Open the left and the right temples.

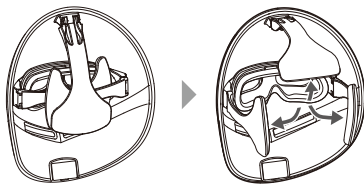HOW TO USE


(Using the Mask and the Neck at the same time)

- ⑦ Wear the Mask on your face.

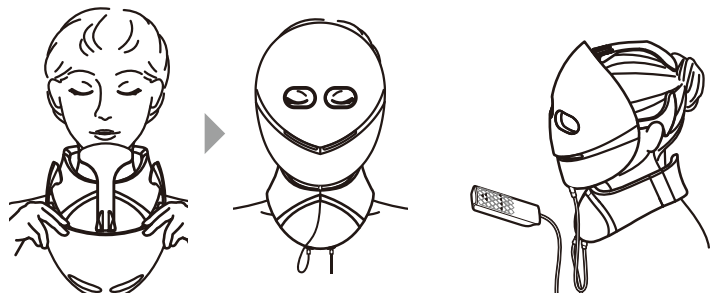
■ Using the Mask Mode 1 or 2

- ⑧ Long press the Button 1 of the remote controller to start the main body in Mask Mode 1.
⑨ Each time you short press the Button 1 of the remote controller, the mode changes from Mask Mode 1 to Mask Mode 2, and to power off. * You can also long press the Button 1 of the remote controller to turn off the power.

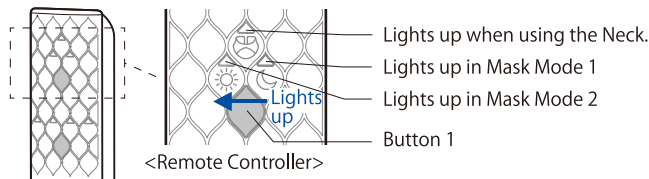■ Using the Mask Mode 3, 4, 5 or 6

- ⑧ Long press the Button 2 of the remote controller to start the main body in Mask Mode 3.
⑨ Each time you short press the Button 2 of the remote controller, the mode changes from Mask Mode 3 to Mask Mode 4, to Mask Mode 5, to Mask Mode 6, and to power off.
* You can also long press the Button 2 of the remote controller to turn off the power.

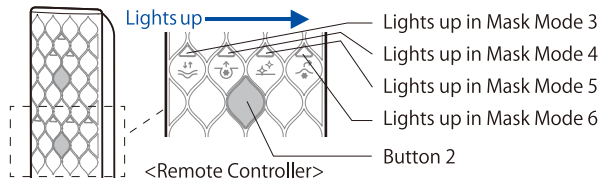● Useful tips

If you long press the Button 2 when using Mask Mode 1 or 2, you can switch the mode to Mask Mode 3.


Likewise, if you long press the Button 1 when using Mask Mode 3, 4, 5 or 6, you can switch the mode to Mask Mode 1.

HOW TO USE (Using the Mask only)


Please make sure this product is fully charged before use. Do not wear your glasses when you use this product.

* You cannot use commercially available USB cables. Use the USB cable provided with the product.


- ① Plug the USB cable (long) into the USB port of the remote controller.

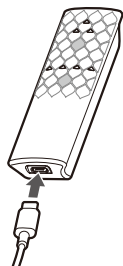
- ② Plug the USB cable (long) into the USB port of the Mask.

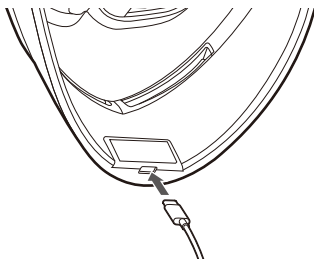- ③ Lift up the head stopper. Open the left and the right temples.

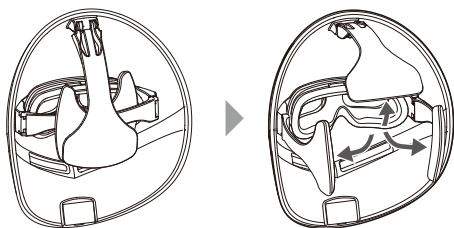
- ④ Wear the Mask on your face.

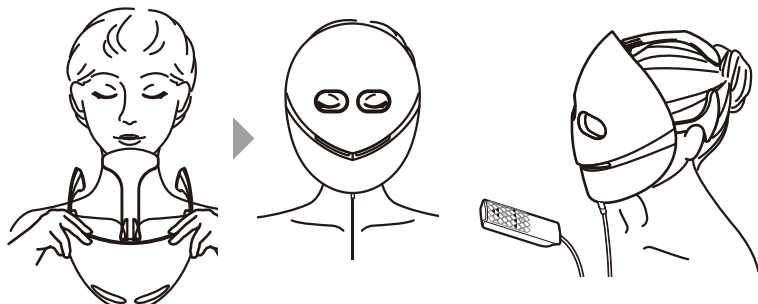HOW TO USE (Using the Mask only)


■ Using the Mask Mode 1 or 2

- ⑤ Long press the Button 1 of the remote controller to start the main body in Mask Mode 1.
- ⑥ Each time you short press the Button 1 of the remote controller, the mode changes from Mask Mode 1 to Mask Mode 2, and to power off. * You can also long press the Button 1 of the remote controller to turn off the power.

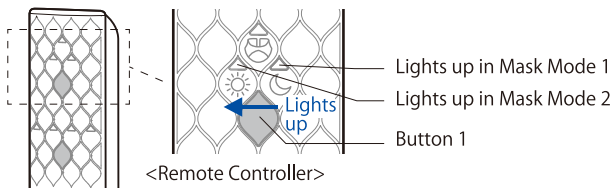■ Using the Mask Mode 3, 4, 5 or 6

- ⑤ Long press the Button 2 of the remote controller to start the main body in Mask Mode 3.
- ⑥ Each time you short press the Button 2 of the remote controller, the mode changes from Mask Mode 3 to Mask Mode 4, to Mask Mode 5, to Mask Mode 6, and to power off. * You can also long press the Button 2 of the remote controller to turn off the power.

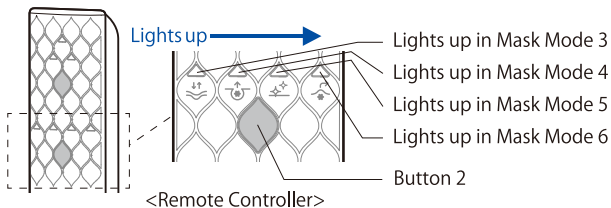● Useful tips


If you long press the Button 2 when using Mask Mode 1 or 2, you can switch the mode to Mask Mode 3.
Likewise, if you long press the Button 1 when using Mask Mode 3, 4, 5 or 6, you can switch the mode to Mask Mode 1.

HOW TO USE (Using the Neck only)


Please make sure this product is fully charged before use.

* You cannot use commercially available USB cables. Use the USB cable provided with the product.

- ① Plug the USB cable (long) into the USB port of the remote controller

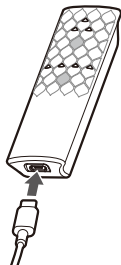
- ② Plug the USB cable (long) into the USB port of the Neck.

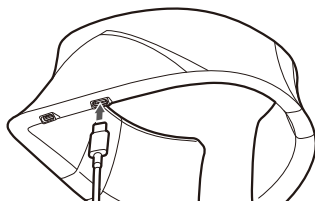
* You can plug the cable into either of the two USB ports.

- ③ Open the Neck with your hands to put it around your neck.

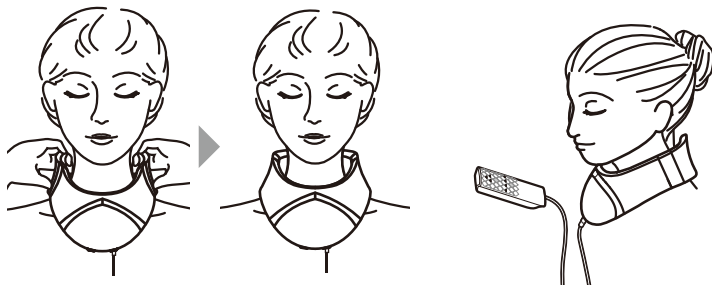- ④ Long press the Button 1 of the remote controller to start the main body in Neck Mode.

- ⑤ Short press the Button 1 of the remote controller to turn off the power.

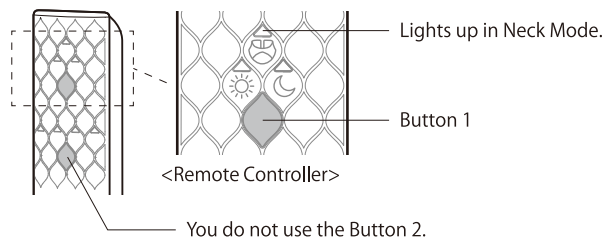
Question	Answer
What is an LED?	LED stands for light-emitting diode, which is a type of semiconductor device that emits light when current flows through it. This product features a linear arrangement of LEDs of different colors that radiate light.
What should I do if I have a skin issue?	Discontinue use of the product until your skin condition heals, and consult a doctor. Skin condition also depends on body condition.
How should I take care of the product after use?	Wipe the product clean using a soft cloth dampened with water. Wiping with a hard cloth may leave scratches on the surface or cause a decline in performance. If the eye pad becomes heavily soiled, you can detach the eye pad to wash it. If you rinse the eye pad, make sure it is completely dry before reattaching it to the main body. * Do not use alkaline agents, alcohol, thinner, benzene or other solvents, polish removers, polishing powder, hot water or chemicals as these can cause deformation or failure.
Can I use the product outside of Japan?	You can use the product outside Japan because it is rechargeable and because the AC adapter can work at voltages of between 100 V and 240 V. Check the voltage and the shape of power outlets in the country where you intend to use the product, and purchase a converter plug if needed.
Can I use commercially available USB cables?	You can use only the USB cable provided with the product. If you have lost or cut the cable, contact your dealer or our customer support center.
Can I replace the battery by myself?	You cannot replace the battery yourself as this may cause failure. Contact your dealer or our customer support center.

<How to detach the eye pad>

Hold the eye pad and pull it.

To reattach the eye pad, press it into the main body.

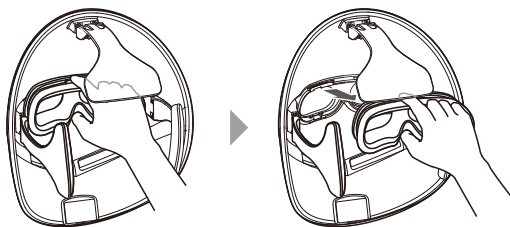⚠ DANGER

When you attach the eye pad, make sure before use that it is fully and firmly attached to the Mask.

TROUBLESHOOTING

Problem	Steps
It does not turn on when I press the button.	<ul style="list-style-type: none">• The remaining battery life may be low. Recharge the battery by following the steps set out in the instruction manual.• Check to make sure the remote controller and the main body are correctly connected.• Check to make sure the dedicated USB cable is used.
When I charge the product, the indicator does not light up.	<ul style="list-style-type: none">• Check to make sure the remote controller and the dedicated AC adapter are correctly connected.• Check to make sure the dedicated AC adapter is correctly plugged into a power outlet.• Check to make sure you are charging using the dedicated AC adapter and the dedicated USB cable.
The product turns off soon after being switched on, even when charged to full capacity.	<ul style="list-style-type: none">• The product may not have been charged correctly. Check to make sure you are charging using the dedicated AC adapter and the dedicated USB cable.• The battery may have degraded. Contact your dealer or our customer support center.

If you cannot solve the problem by following the steps above, contact your dealer or our customer support center.

PRODUCT SPECIFICATIONS

Product name	CHOUOHC THE RADIANCE KIT (Mask & Neck)
Model number	CRMN-(Color Number)-01

Product name	CHOUOHC THE RADIANCE MASK
Model number	CRM-(Color Number)-01

Product name	CHOUOHC THE RADIANCE NECK
Model number	CRN-(Color Number)-01

Main Body

Battery	Rechargeable
Built-in battery	Lithium-ion battery 3.6 V 5,000 mAh
Charging time	about 3 hours
Number of LEDs	Mask : 1026 LEDs Neck : 502 LEDs
Weight	Mask : approx. 360g Neck : approx. 190g Remote controller : approx. 140g
Size	Mask : W201 × H254 × D209mm Neck : W174 × H90 × D207mm Remote controller : W37 × H130 × D35mm
Intended use	Household, indoor

AC Adapter

Rated input voltage	AC100-240V
Rated frequency	50/60Hz
Rated output	5V 2A
Weight	approx. 40g
Size	W36.6 × H23.8 × D40.5mm

Country of Manufacture

Mask : Japan

Neck : Japan

Remote controller: China

RECYCLING OF LITHIUM-ION BATTERY

Make sure to take out the lithium ion battery when disposing the device.

Make sure to consume the power in the lithium ion battery before taking out the lithium ion battery.

Do not discard the lithium ion battery. Please follow your local government rules to dispose in the specified battery recycling place.

Li-ion

(Battery removal procedure: Unclench the shell of the device with a tool, then take out the battery. Be careful to avoid the battery's positive and negative short-circuiting)

ASSURANCE PROVISION

This warranty guarantees free repair subject to the conditions included in the instruction manual. If the product is damaged within the warranty period from the date of purchase, please request a repair from your dealer or our customer support center with a copy of the warranty attached.

* The warranty is not valid if the serial number, date of purchase, name of dealer and customer name, address and phone number are not provided. Please make sure in advance that you have all the required information.

* If the serial number sticker on the warranty is removed, the product is not covered by the warranty.

◆ Repair costs will be charged in the following cases even if the product is within the warranty period:

- A. The serial number sticker on the warranty has been removed
- B. The serial number, date of purchase, name of dealer or customer name, address and phone number are not provided in the warranty
- C. Damage or failure has occurred due to improper handling (e.g. dropping and/or submerging in water) during transportation or moving by the user
- D. Damage or failure has occurred due to fire, natural disaster or abnormal voltage
- E. Failure has occurred due to the connection of other equipment to this product
- F. Failure has occurred due to the use of an AC adapter or USB cables other than the dedicated items provided with the product
- G. The product has been used incorrectly or repair, adjustment or modification has been performed by anyone other than our company
- H. The product is used for non-household use (e.g. commercial use)
- I. The warranty is not submitted
- J. Damage has occurred due to the use of any beauty essence

◆ This warranty is valid only in Japan.

* In cases where the product is purchased outside of Japan, the product is covered by warranty only where an authorized distributor certificate has been issued by the manufacturer.

◆ Shipping fee

- a. The company will bear delivery and return costs during the warranty period.
- b. Delivery and return costs will be borne by the customer after the warranty period.

This warranty will not be reissued. Please store it in safe place.

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